


PLUS

WWW.ROTAID.COM

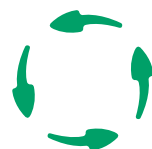

ROTAID™
AED CABINETS

ROTAID PLUS

PRODUCT SHEET

The Rotaïd Plus Heat is an innovative AED wall cabinet that ensures a recognizable and sustainable storage of your AED device. The characteristic round shape and green colour, make it easy to be found and recognized and comes with a high quality appearance. The large AED logo, positioned in the transparent centre of the cover, leaves no room for interpretation. This is a life saving device.



The cabinet can be opened by an intuitive twist that removes the whole cover. The responder therefore has full access to the AED that hangs via a bracket for a quick and easy grab. The cover is attached to the back plate via a belt and needs no further attention from the responder. A stainless steel hook, attached to the bottom of the cabinet, can hold a seal, tag or other appliance.

The cover, equipped with recesses for an optimal grip and directive arrows, is made from 100% Polycarbonate. With a thickness of 5 millimetres it offers a solid protection and ensures a high grade of water- and dust resistance. Furthermore, the cover is fitted with UV filters, to protect from direct and indirect sunrays.




The Plus can optionally be equipped with a Visual and acoustic alarm and AED standard (to use with AED without handle)




By modern production techniques and the use of recycled plastics, the Plus meets the highest requirements when it comes to sustainability.

ROTAID PLUS

TECHNICAL DETAILS

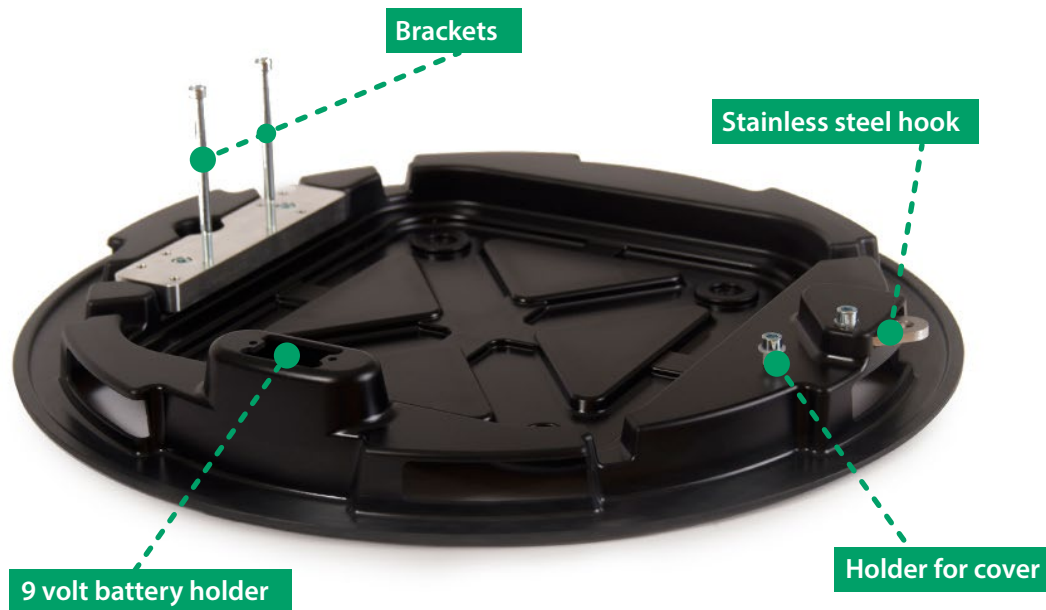

The back plate is made from 100% ABS plastics and offers space and protection to practically all types of AED devices on the market.

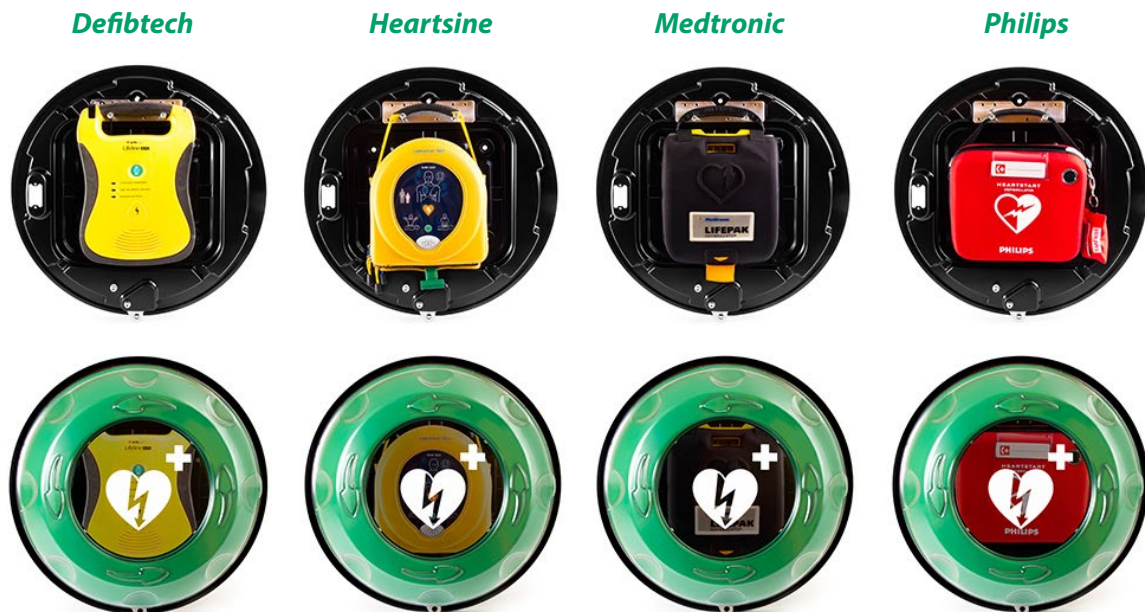

Also compatible with *Cardiac Science, CU Medical, Primedic, Schiller and Zoll.*

Dimensions

Length **47 cm** Width **47 cm** Depth **18 cm** Weight **1,7 kg**

Length **18,5"** Width **18,5"** Depth **7,1"** Weight **3,75 lbs**

Packed individually by box, fastening materials included.